


# Pb Shipper Now Available for Cardiac Size Flood Source

Save time and money with our “hassle-free” Pb Shipper. With a lead-lined sleeve covering all 6 sides, our new Pb Shipper is designed to qualify as a UN2910 “Excepted Package” when the surface reading for your Co-57 flood source is  $\leq 0.5\text{mRem/hr}$ .

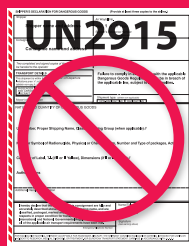

Do you find yourself holding on to your depleted flood sources longer than you want?

Pb Shipper is the convenient and efficient alternative for returning and receiving your flood sources.

 **Eckert & Ziegler**  
Isotope Products

\* The trained shipper who is offering the package for shipment is responsible for determining the proper shipping name of the package. If the packaged source exceeds  $0.5\text{mRem/hr}$  at the surface of the package, you will need to ship the source as a Type A UN2915 package – shipping papers and emergency information documentation will need to be included with the UN2915 shipment.



## ProKem Technology™ = Better Quality Control

ProKem Technology begins with our proprietary chemistry which results in a unique purification and production process that produces cleaner Cobalt for a more uniform source. This yields superior camera calibration, less variation and, ultimately, better quality control.



**Eckert & Ziegler**

Isotope Products

**Eckert & Ziegler Isotope Products**  
24937 Avenue Tibbitts  
Valencia, CA 91355

### USA

Phone: +1 (661) 309-1010  
Fax: +1 (661) 257-8303  
Toll Free: +1 (866) 476-9767

nucmedsales@ezag.com  
www.ezag.com

**Eckert & Ziegler Isotope Products GmbH**  
Niederlassung Braunschweig  
Harxbütteler Str. 3  
38110 Braunschweig

### Germany

Phone: +49 5307 20165-60  
Fax: +49 5307 20165-666

ipe@ezag.de  
www.ezag.com

**Eckert & Ziegler Isotope Products Sarl**  
12 Avenue des Tropiques  
Hightec Sud - Bâtiment B  
91955 Courtaboeuf Cedex

### France

Phone: +33 1 64 86 22 22  
Fax: +33 1 69 86 10 65

jeanrene.caramelle@ezag-france.fr  
www.ezag.com